



## USER MANUAL



Use it at home.



Use it in a hotel whilst on holiday.



Take it on camping trips.



Use it on cruise ships.



Take it with you when visiting friends.

1

Lift the toilet seat and slide the pipe attached to the spray arm through the seat fixing posts.



2

Gently press down on the spray arm so the soft suction cups fit firmly to the back of the basin.



3

Connect the box and douche arm firstly by attaching the small pipe into the Porta-Bidet unit.



4

Next, attach the short pipe from the Porta-Bidet to the long pipe attached to the spray arm.




5

Connect the control button to the Porta-Bidet by inserting the 3.5mm jack to the underside of the bidet.



6

Add a litre of water to the bidet's reservoir at your preferred temperature using the filler jug provided.

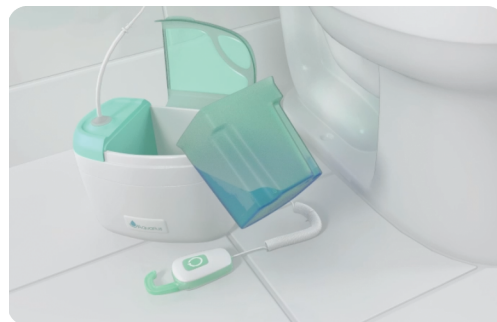

7

Press the button to activate the pulsating flow.



## SPRAY SETTING

Setting the spray arm pressure for your own comfort is easy.

Simply set up the Aquarius Porta-Bidet in a normal manner, then hold down the button on the control unit and the Porta-Bidet will reset itself starting with a light trickle. The water pressure will continue to rise. Once your comfort level is reached, release the button on the control unit.

The Aquarius Porta-Bidet will now memorize this setting for your future use.



PLEASE KEEP THESE INSTRUCTIONS!

## Charging

On first receiving your Aquarius Porta-Bidet, place on charge for at least 12 hours. This will condition the battery.

Place the charger connector into the charging point on the underside of the Aquarius Porta-Bidet, marked *Charger*.

Please note the unit will not operate when charging.



### WARRANTY DETAILS

The Aquarius Porta-Bidet has a 12 month return to base warranty.

Please see website for details:

[www.aquariushygiene.com](http://www.aquariushygiene.com)

If you have questions, please contact our customer support line:

Tel no. 0800 064 4454 (option 4)

Or email us at:

[support@aquariushygiene.com](mailto:support@aquariushygiene.com)